



Medical Therapy Jaundice PL-L – the compact way to eliminate the need for blood transfusions

Blau (/52) – PL-L, PL-S

By emitting light almost entirely within the 400 to 500 nm bandwidth these Medical Therapy Jaundice PL-L lamps have no radiation from the short wave UVB waveband. They are therefore ideal for treating new born babies suffering from hyperbilirubinemia (neonatal jaundice) and Crigler-Najjar Syndrome (CNS). Moreover, the bandwidth of these lamps peak at the most effective treatment wavelength of 450 nm. This highly efficacious phototherapy treatment has eliminated the need for blood transfusions in almost all jaundiced infants. In addition, with the compact format of this PL-L lamp equipment manufacturers have more design freedom in developing their solutions. By emitting a full spectrum of a high color temperature, the PL-L /953 lamps are ideal for treating SAD.

Vorteile

- Optimales Spektrum für verschiedene Anwendungen in der Phototherapie

Merkmale

- Emissionsspitze bei 450nm

Anwendung

- Phototherapie
- Landwirtschaft

Blau (/52) - PL-L, PL-S

Hinweise


- Es ist sehr unwahrscheinlich, dass der Glasbruch einer Lampe negative Auswirkungen auf Ihre Gesundheit hat. Wenn es zu einem Glasbruch kommt, lüften Sie den Raum 30 Minuten lang und entfernen Sie die Splitter möglichst mit Handschuhen. Legen Sie die Splitter in einen verschließbaren Plastikbeutel und geben Sie ihn beim Werkstoffhof zum Recycling ab. Benutzen Sie keinen Staubsauger.

Versions



XPPR XUMPLL 2G11

Abmessungsskizzen

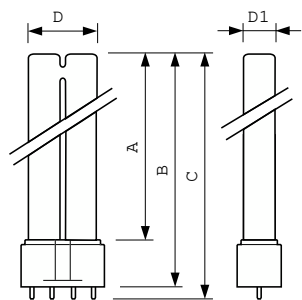

Product	D1 (max)	D (max)	A (max)	B (max)	C (max)
PL-L 18W/52/4P 1CT/25	18 mm	39 mm	188,2 mm	214 mm	220,6 mm

